


Navigating the SC 12-Inch Sash Menu

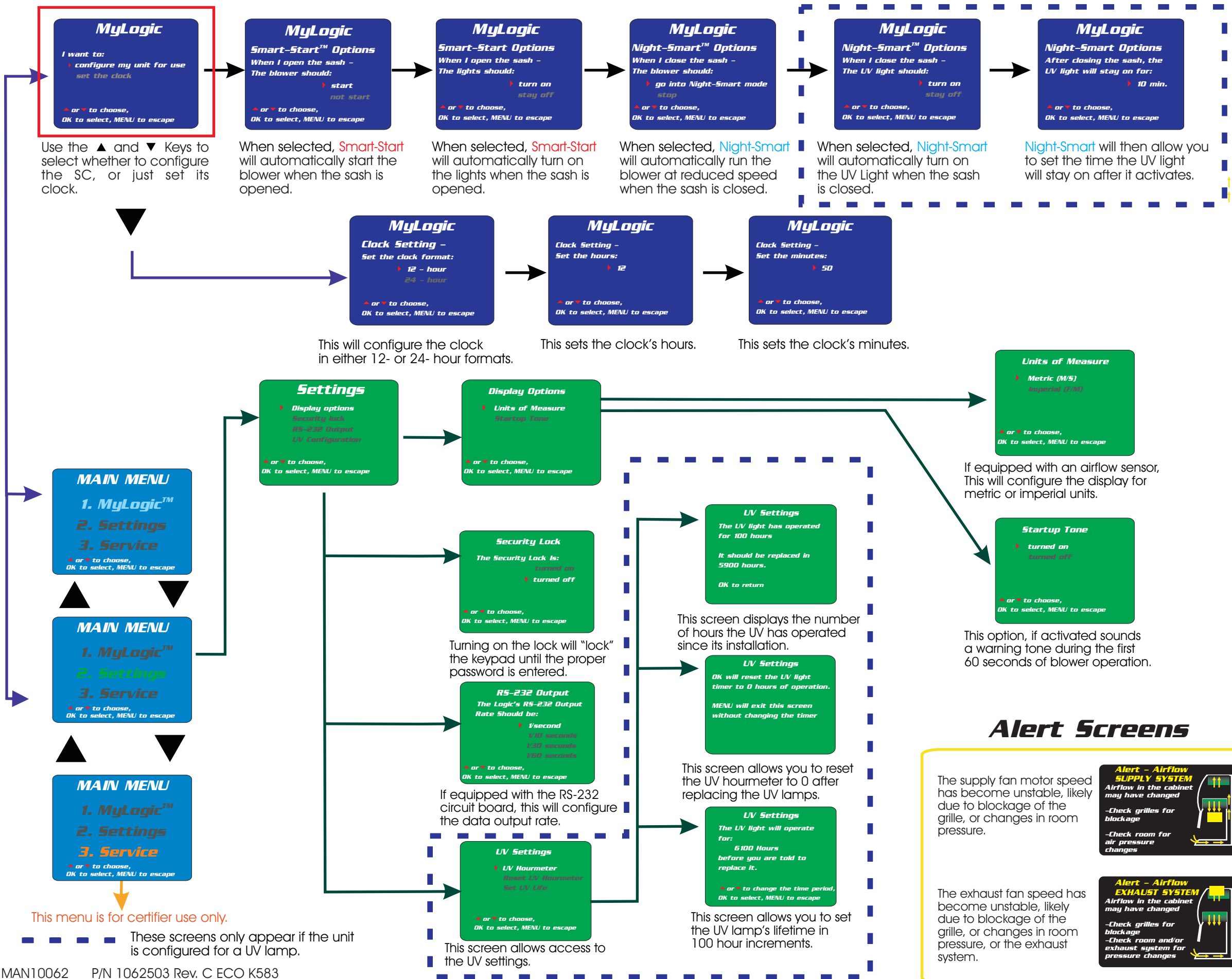

Alarm Screens

Alarm - Sash Height
The sash is open too far for safe operation.
-Lower it
-Press Mute/OK to silence the alarm for 5 minutes

Alarm - Airflow
The Inflow readings are out of range-
-Have a certifier recalibrate the airflow sensor

Alarm - Exhaust Fan
The exhaust fan has failed
-Have a Certifier examine the cabinet

Alarm - Exhaust System
There is insufficient exhaust flow -
ACTIVE PROTECTION ENGAGED

Alarm - System
There is a loss of communication in the control system-
-Have a certifier examine the cabinet

Alert Screens

Alert - Airflow SUPPLY SYSTEM
Airflow in the cabinet may have changed
-Check grilles for blockage
-Check room for air pressure changes

Alert - Airflow EXHAUST SYSTEM
Airflow in the cabinet may have changed
-Check grilles for blockage
-Check room and/or exhaust system for pressure changes

Warm up Screen

PLEASE WAIT!

This screen is displayed in ALL configurations during the first 60 seconds of operation.

Navigating the SC 12-Inch Sash Menu

Diagnostics

